

Device Summary

LUCAS® 3 Chest Compression System


Incident ID:
 Patient ID:
 Patient Name:
 Power On: 5/24/2018 10:39:40 AM
 Recording Duration: 00:05:10
 Compressions Duration: 00:04:59

Device Type: LUCAS® 3
 Device ID: LUCAS3-35160044-ANI
 Serial Number: 35160044
 Software Version: LUCAS3 v3.1.4.26


*Times have been adjusted by the system.

LUCAS Data

LUCAS compression ratio

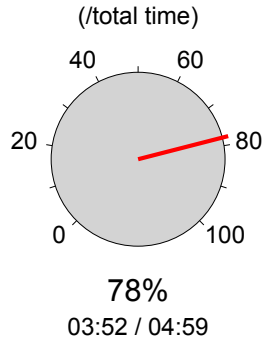

Compression rate

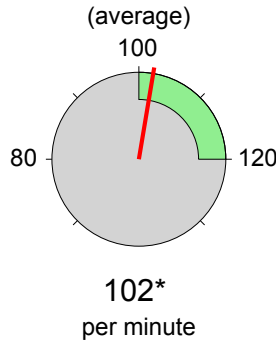

Compression count: 399
 Pauses over 10 sec: 1
 Longest pause: 14 sec

*Rate changed by user during operation

LUCAS Device Configuration

Compression

Rate at Power ON
 102 111 120

Rate change: ON

Allow rate to alter between
 102 111 120

Depth (nominal patient)

Chest height > 7.3 in: 2.1 ± 0.1 in
 Chest height 6.7 to 7.3 in: 1.5 to 2.1 ± 0.1 in

Ventilation

In 30:2 mode

Ventilation pause: 3.2 sec

Compression/ventilation ratio

30:2 50:2

In continuous mode

Audible alerts: OFF

Alerts/minute: 10

Ventilation pause: 0.3 sec

Suction Cup Start Position

Use AutoFit
 Use QuickFit (1.2 in)
 Use Manual

Pressure Pad Release

Return up in ADJUST mode
 During ventilation pauses (0.4 in)
 In PAUSE mode (0.4 in)
 In ACTIVE mode (0.4 in)

Audible Timer

OFF
 CPR timer
 Continuous timer

Data Transmission

Auto transmit at charging

Device Summary


LUCAS® 3 Chest Compression System

Power On: 5/24/2018 10:39:40 AM

Serial Number: 35160044

*Times have been adjusted by the system.

LUCAS Time Graph

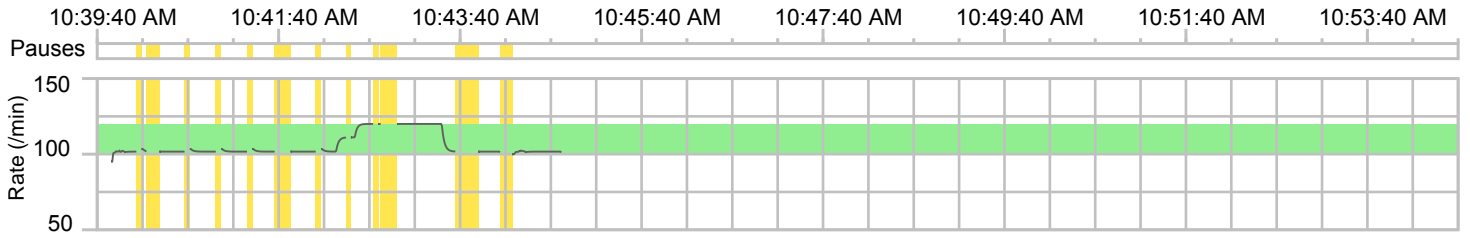

Event log LUCAS® 3 Chest Compression System

Incident ID:	Device Type: LUCAS® 3
Patient ID:	Device ID: LUCAS3-35160044-ANI
Patient Name:	Serial Number: 35160044
Power On: 5/24/2018 10:39:40 AM	Software Version: LUCAS3 v3.1.4.26
Recording Duration: 00:05:10	
Compressions Duration: 00:04:59	

Elapsed Time	Real Time	Description
Thursday, May 24, 2018		
00:00:00	10:39:40 AM	LUCAS 3 Power On
00:00:00	10:39:40 AM	Adjust Mode
00:00:05	10:39:45 AM	Pause Mode
00:00:09	10:39:49 AM	30:2 Mode
00:00:09	10:39:49 AM	First LUCAS Compression
00:00:33	10:40:13 AM	Pause Mode
00:00:35	10:40:15 AM	Adjust Mode
00:00:39	10:40:19 AM	Pause Mode
00:00:40	10:40:20 AM	30:2 Mode
00:01:58	10:41:38 AM	Pause Mode
00:02:07	10:41:47 AM	30:2 Mode
00:02:38	10:42:18 AM	Rate Change to 111
00:02:50	10:42:30 AM	Rate Change to 120
00:03:08	10:42:48 AM	Pause Mode
00:03:18	10:42:58 AM	Continuous Mode
00:03:48	10:43:28 AM	Rate Change to 102
00:03:58	10:43:38 AM	Pause Mode
00:04:11	10:43:51 AM	Continuous Mode
00:04:27	10:44:07 AM	Pause Mode
00:04:29	10:44:09 AM	LUCAS 3 Power Off
00:04:32	10:44:12 AM	LUCAS 3 Power On
00:04:32	10:44:12 AM	Continuous Mode
00:04:32	10:44:12 AM	Battery Replaced 34170305489->34170707468
00:04:56	10:44:36 AM	External Power Connected
00:04:58	10:44:38 AM	External Power Disconnected
00:05:00	10:44:40 AM	External Power Connected
00:05:07	10:44:47 AM	Pause Mode
00:05:07	10:44:47 AM	Last LUCAS Compression
00:05:10	10:44:49 AM	LUCAS 3 Power Off